

TAURUS 






Continental delivery type



International delivery type

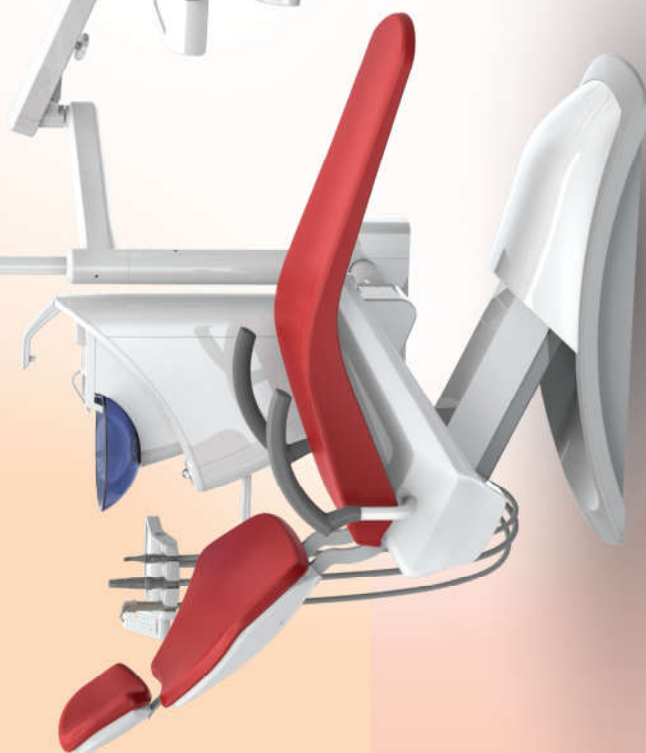

TAURUS 

Minimal design

Chair mounted dental unit
Ergonomic design offers a space saving and very efficient working environment



simple
design

Basic configuration

Selected essential configurations are equipped
Many parts are removable for disinfection
Light sensor cup filler




essential
function

Affordable pricing

Reasonably priced, equipping only necessary configurations
Various parts and instruments are offered economically



reasonable
price

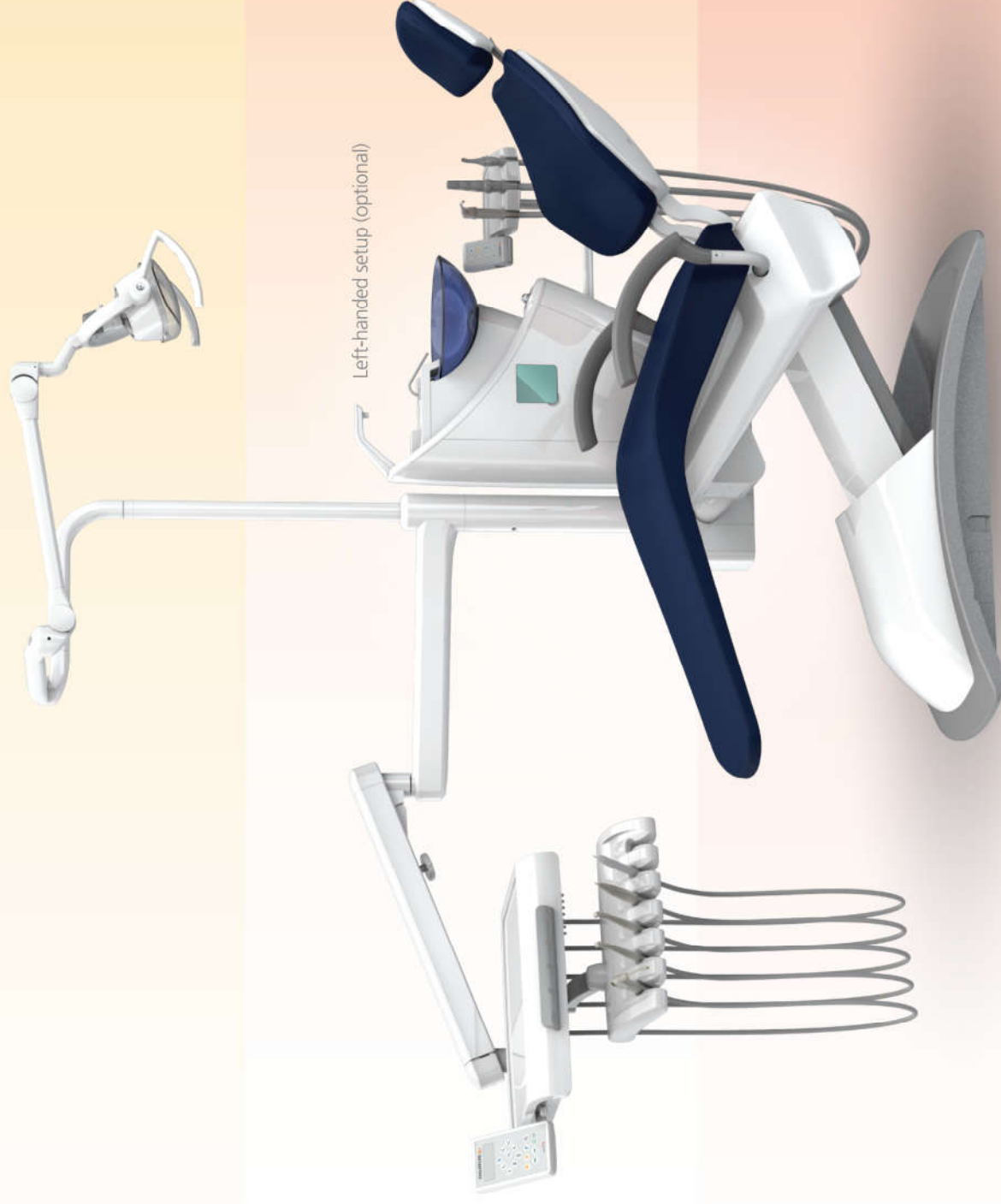

Left-handed setup (optional)

Spittoon

The spittoon bowl is easily removable without the help of any instrument and easy to clean and wash



Twin-articulating headrest

Various angles are enabled on the headrest, providing optimum positions for both adults and children (optional)



Operating light

The light reaches to optimal positions for both maxillary and mandibular treatments. (C-200 LED operating light is optional)




Rotating Armrest

The armrest ensures zero obstruction to patient entry and exit ; also provides comfort for the patient's arm


Safety Auto System

Safety feature of a backrest and assistant table ensures retraction if contact is made with any object

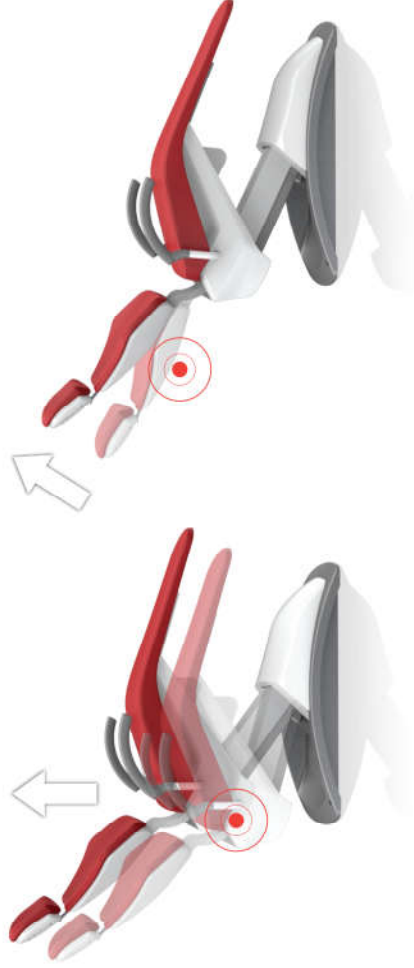

Rotating chair-mounted cuspidor unit


The compact swiveling body of the device allows the assistant to have better access to the patient

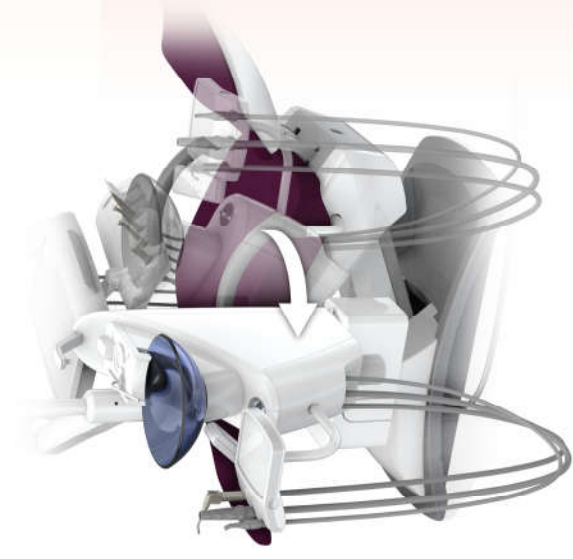

Seat Upholstery




Red

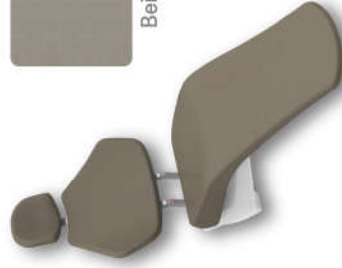

Beige



Bronze

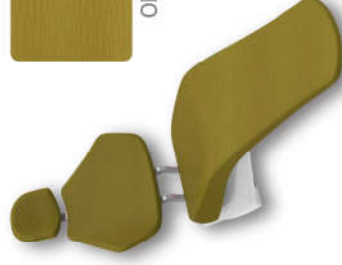

Olive



Blue



Navy



Purple



Black



Orange



Mango



Golden brown



Mint green



Dark Blue



Purple pink



Walnut

Specification

Option




EMS LED Scaler



MC3



LED Light

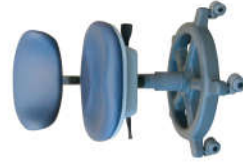

Dentist Stool



Assistant stool



Optima MX2 INT

Standard



Saliva Ejector



Vacuum Syringe



CL-2000 Light



Foot Controller

Rated voltage : AC 220V, 60Hz (110V, 50Hz available)

Power consumption : 1.0kVA

Seat movement powered by hydraulic motor

Chair elevation range : Max. 695±20mm, Min. 400±20mm

Backrest inclination : 0±2° ~ 70±2°

Seat length : 1820mm

Weight : 165±5kg

Chair auto position program : 2 positions

We reserve the right to modify the design or specifications with prior notice to allow for the improvement.



Shinhung. Co., Ltd.

450, Cheongpa-ro, Jung-gu, Seoul, Korea

Tel : 82-2-6366-2175 / Fax : 82-2-6366-2131

<http://www.shinhung.co.kr>

E-mail : export@shinhung.co.kr

We reserve the rights to modify the design or specification without prior notice to allow for the improvement.

CC/1511/2000 Rev. 0