



MINIMIZE DOUBT  
MAXIMIZE SATISFACTION



## HEAD LIGHT & CAMERA

### SHADOWLESS HEAD LIGHT

- + NO SHADOW FROM MEDICAL INSTRUMENTS
- + BRIGHTER THAN THE CONVENTIONAL HEAD LIGHT
- + SEMI-PERMANENT LIFE OF LIGHT SOURCES
- + LIGHT WEIGHT (ABOUT 100GRAM) & WIRELESS
- + EASY ON & OFF
- + INSTALLABLE CUSTOMIZED LOUPE OR SURGICAL TELESCOPE
- + HIGH QUALITY LI-ION RECHARGEABLE AND BATTERIES



DR.KIM | HEAD LIGHT  
DKH-50


LIGHT ON



LIGHT OFF



CONVENTIONAL HEAD LIGHT

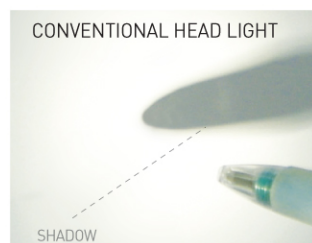

SHADOW

DR.KIM HEAD LIGHT



SHADOW

2 Light beams from 2 LED light sources converge in a point 30~40cm from the eyes. This head light reduces the shadow of medical instruments

## HEAD LIGHT

• **DKH 50** • REMARKS : CHARGER & 2 RECHARGEABLE BATTERIES INCLUDED



### The world first ergonomic C-BAND Type

- + Illumination: Variable (0~70,000lux)
- + Operation continuity : 2.5~12hours
- + Weight(Exception rechargeable battery: 51g): 104g
- + Affordable head length : 160-210(mm)
- + The angle of illumination from horizontal to downward : 0-41 degrees
- + Included 1 charger and 2 rechargeable batteries
- + Charging time : about 2 hours
- + Use Dr. Kim's specific charger and batteries for safety
- + Installable customized loupes or surgical telescope

## CAMERA

• **D-CAM2** • ATTACHABLE TO HEAD LIGHT (DKH-50)



### Auto Focus Camera

- + Sensor : 5M Pixel
- + Interface : USB 2.0
- + Weight : 13g
- + Resolution : Full HD, HD
- + Auto Focus / Manual Focus
- + 2X Digital Zoom
- + Attachable to AKH-50
- + Real-time Display & Recording
- + Operating in Window OS



Suite 502, 15 Orion Rd, Lane Cove, NSW 2066

T 02 8084 2900 F 02 8084 2700

[www.minimaximplant.com.au](http://www.minimaximplant.com.au)

Email [info@minimaximplant.com.au](mailto:info@minimaximplant.com.au)

1 year

WARRANTY  
FROM SHIPPING DATE



RETURN  
GUARANTEE